
currentbody.com

Version 1 - August 2023

CURRENTBODY|skin Laser Hair Removal Device is manufactured for
CurrentBody Ltd, Q17, Crossley Road, Stockport, Greater Manchester, SK4 5BB, United Kingdom
EU Authorised Representative: Certification Company, Veluwezoom 42, 1327 AH, ALMERE, The Netherlands

CURRENTBODY|skin™

LASER HAIR REMOVAL DEVICE

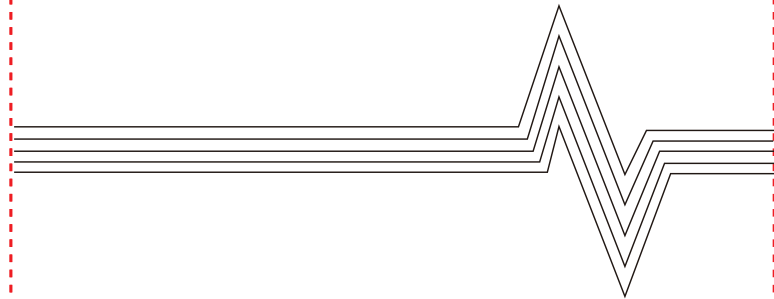USER GUIDE AND MANUAL

Includes important safety information.

Read all instructions prior to use.

PLEASE KEEP THESE INSTRUCTIONS SAFE

CONTENTS

CurrentBody Skin™ Laser Hair Removal Device	02
Device Overview	03
What is in the box?	04
Important Safety Information	05
Suitability Chart	06
Warnings	07
Cautions	10
How to Use	11
Maintenance	15
Warranty	16
Troubleshooting	17
Symbols explained	18
Technical Specification	19
French	20
Italian	36
German	52
Spanish	68

Welcome to CurrentBody Skin™ Laser Hair Removal Device

Established in 2009, CurrentBody is the world's only online retailer dedicated to home-use beauty devices. As the first company to recognise devices as one of the fastest growing trends in the world of beauty, we have unrivalled expertise in selecting at-home devices that are science-backed and clinically proven, offering our customers the ultimate choice of devices that treat a wide range of concerns, including anti-ageing, cleansing and hair removal.

Underpinning our range of innovative beauty technology is expert education to help our customers choose the right device for them and excellent customer service.

Our highly trained beauty device experts are on hand to offer impartial advice and answer your questions in person, as well as how-to videos and top tips available throughout our website and social media.

DEVICE OVERVIEW

Using professional technology found in clinics and salons around the world, the CurrentBody Skin™ Laser Hair Removal device is one of the most powerful hair reduction systems available in the home-use market. Using laser technology, you can experience up to 98% hair growth reduction in just 12 weeks.

Suitable for use on the face, underarms, bikini line, arms, back and legs, you can treat the whole body in as little as 10 minutes. The concentrated beams of laser energy target the pigment (melanin) in the hair follicles, inhibiting their ability to regrow hair.

Week 1-3 Use 3x per week

Week 4-12 Use 1x per week

Week 12+ Use 1x per month

It is intended for use only on light to medium skin tones (see the skin tone chart) and naturally brown or black hair. It will not work on white, grey, blonde or red hair.

WHAT IS IN THE BOX

CURRENTBODY Skin™
Laser Hair Removal
Device


Power adapter

Power cord

Important features:

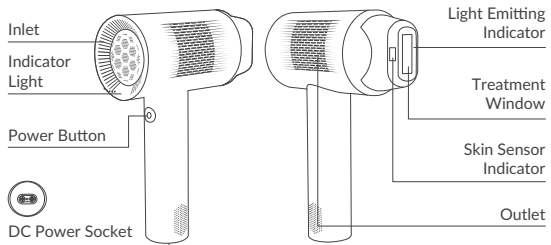IMPORTANT SAFETY INFORMATION

- DO NOT operate the device where it could come in contact with water.
- DO NOT use the device with a power adapter that is damaged or faulty. Do not harshly bend or abuse the power cord.
- DO NOT attempt to take the device apart. The device is not serviceable by the user. This will void the warranty.
- STOP using the device immediately if you think it is defective.

SUITABILITY CHART

To check if your skin tone is suitable for use on the laser hair removal device, consult the suitability chart below. You may be injured (burnt, discolored, or scarred) if your skin tone is not listed as an intended skin tone for the treatment.

This device is designed to work on naturally light brown to black hair. Lighter hair colours such as white, blonde, red or grey will not absorb enough of the laser's energy to disable the hair follicle, meaning the treatment will not work.

Suitable Hair Colour

Suitable Skin Tone

For Use On

Legs, arms, bikini line, underarms, feet, hands, chest, back, stomach, women's facial hair (cheek line downward).

Not suitable for use on a man's face, jaw, neck and genitals.

Customer Care

Please contact CurrentBody via www.currentbody.com if you have any questions or need any advice.

WARNINGS

- This device is designed for use on adults over 18 years of age and not recommended for use on infants or children.
- Please consult your doctor prior to use if you are, or suspect you are, pregnant or are breastfeeding.
- Please consult your doctor prior to use if you have or suspect you have cancer.
- Do not use the laser on a man's face, jaw or neck. The hair in a man's beard area is too dense, and using the laser on these areas may result in skin injury.
- Do not use if your skin is red or swollen from a previous treatment.
- Do not look directly into the laser. An eye injury may result from misuse.
- Do not use the device to treat skin that is painful to touch, raised (higher than surrounding skin), crusting, or bleeding. This may be a sign of a more serious condition and should be checked by your doctor.
- The device should not be used if you suffer from diabetes mellitus, coagulation disorders, or malignant skin lesions.
- This device should not be used on areas such as birthmarks, moles, tattoos and warts, and should not be used after exposure to the sun.
- If the Suitability Chart indicates that your skin tone is not suitable for the device, you should not use it.
- It is recommended that certain medications, such as blood thinners or topical medications such as Retinol, be avoided prior to laser hair removal.
- Do not use medications or topical treatments which cause photosensitivity. Please check with your doctor before use.

- It is not recommended to use the device directly on the ears, nipples, genitals, or around the anus. In these areas, the skin may be more sensitive, the skin color may be darker, and/or the hair density may be greater. There is the possibility of serious skin injuries resulting from the use of the device there.
- It is not recommended to use the device if you have skin conditions such as tinea, eczema, herpes, or other dermatitis.
- Do not leave the device still for more than 8 seconds on the same area of the body.
- No modification of this equipment is allowed.
- Before disposal of this device, please review your local laws and requirements surrounding disposal of electrical items. The preferred method of disposal is through recycling.
- The device should not be used near bathtubs, showers, basins, or other water-containing vessels. Do not use the device if it is wet.
- During treatment, ensure that the air outlet is not blocked.
- Place the device away from heat-generating devices or humidifiers.
- A flammable or explosive environment should not be used with this device.
- Do not use the device if it is damaged, for example if the cable is damaged, if the components are cracked or broken, etc.
- Attempting to open or remove the covers of this product may result in fire, electric shock, or damage to the equipment and invalidate your warranty.

- Please use the enclosed power cord to connect to the standard power outlet in your country. Ensure that the power cord is operating within its rated voltage. The failure to do so may result in a fire or electric shock.
- Pull the plug firmly to disconnect the power cord. There is a risk of fire or electric shock if you tug on the cord.
- Be sure to use the correct voltage. A specific voltage is required for the use of the unit. An incorrect voltage connection may result in fire, electric shock, or damage to the device. Overloading a power circuit can result in a fire or electric shock.

CAUTIONS


- Please avoid touching the plug with wet hands.
- It is recommended that you unplug the unit before cleaning it. There is a possibility of electric shock if the unit is being cleaned while plugged into a power outlet.
- Periodically clean the area around the plug with a dry cloth. The plug may catch fire if it is covered in dust, water, or oil.
- After turning off the power switch, disconnect the power cord from the socket if you intend to leave the unit unused for an extended period.
- Discontinue use if the device overheats or becomes uncomfortably hot.
- Use only the supplied power cord and power adapter. Use of an incorrect power cord could damage the device.
- Do not store the device in direct sunlight or in a location with temperatures over 100 Fahrenheit as this may cause damage to the plastic coating.
- Do not attempt to open or repair the device or power supply.
- Do not drop the device on hard surfaces, as this could damage it. Using a damaged device could result in injury.
- If you have any concerns about the skin you are treating, you should cease use and check with your doctor.
- Contact your doctor if you have any medical questions regarding use of this device.

HOW TO USE

Before the First Treatment

- Do not wax or pluck your hair 6 weeks prior to your treatment, and do not undergo a skin peel during this period.
- A patch test should be performed at least one day prior to treatment (see below).

Before Each Treatment

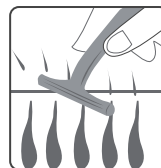a) Ensure skin is clean and dry before use. Remove all cosmetics, lotions, and creams from the treatment area that may affect the laser's performance.

b) Carry out treatment within 4 hours of shaving the area. Removing the surface hair improves comfort and allows the light to better reach the hair in the follicle.

How to Carry Out A Patch Test

Before you start your test, we recommend you check your skin tone against the Skin Tone Suitability chart.

You will also need to shave the area before patch test.

- Press and hold the power button for 2 seconds to turn the device on.
- The device will take a few moments to initialise, then it will enter 'Standby' Mode.
- Press the treatment window against the patch of skin you are wanting to test
- If your skin tone is compatible with the device, the white LED indicator light will illuminate and turn white and the skin tone sensor will start detecting the skin tone. If your skin tone is not compatible, the device will continually 'beep' and the

indicator light of the light-emitting window will flash blue until it leaves the skin.

- If your skin tone is compatible you can now proceed with your treatment.
- Double-press the power button to unlock the laser. The device is in 'Ready' Mode.
- You can now apply the treatment window to the desired area.
- Please carry out one single laser pulse on treatment level 1 (default). If this is tolerable, please move to a new area and try level 2 to establish your ideal comfort level.
- Please increase the level, until you find your desired comfort level.
- Please wait 1 hour after completing your patch test to ensure it is safe to use
- Each time you turn on the device, the skin tone sensor only detects the skin tone once.

Treatment Steps

You will need to shave the area before carrying out a treatment.

- Press and hold the power button for 2 seconds to turn the device on.
- The device will take a few moments to initialise.
- Press the treatment window against the patch of skin you are wanting to treat.
- If your skin tone is compatible with the device, the white LED

indicator light will illuminate and turn white and the skin tone sensor will start detecting the skin tone. If it is not compatible, then it will not turn white.

- If your skin tone is compatible you can now proceed with your treatment.
- Double-press the power button to unlock the laser. The device is in 'Ready' Mode.
- You can now apply the treatment window to the desired area.
- The skin tone sensor indicator will turn white when in contact with the skin. The light emitting indicator will turn from blue to pink/red to indicate that the user's skin matches the condition of hair removal.
- Move the device in a 'stamp and release' motion along the treatment area, moving after each beep and pink/red light indication. Continue with this motion, ensuring there are no overlaps or missed areas. When the treatment window leaves the skin, the emission ceases within 0.1 seconds.
- After the treatment has been completed, hold the power button for two seconds to turn off the device, and then unplug it. Do not unplug the power plug directly without first turning off the device to avoid damage. If the device remains in the 'standby' mode for five minutes without any activity, it will automatically turn off.

Between Treatments

It is recommended that you shave rather than wax or pluck your hair in between treatment sessions.

Continue To Treat

It is recommended that the device is used:

Week 1-3 Use 3x per week.

Week 4-12 Use 1x per week

Week 12+ Use 1x per month

Treatment Maintenance

Once you have achieved your desired results, continue to treat monthly after 12 weeks to achieve and maintain your desired results.

MAINTENANCE

We recommend keeping the components of the device clean, in particular the treatment window and the skin sensor. Please follow the instructions and warnings below when cleaning your device.

Instructions For Cleaning

- The first step is to ensure that the device is turned off. Afterwards, disconnect the power adapter from the device and the power outlet.
- You can clean the device and the power adapter's exterior with a damp cloth. Make sure they are completely dry before applying your next treatment.
- Be careful not to damage the treatment window when removing dirt or debris.
- Wipe the skin sensor gently with a dry, lint-free cloth.

WARRANTY

Under normal use this product is guaranteed for 2 years after its original purchase date against manufacturing defects in materials or workmanship. This warranty extends only to the original purchaser. If the product should become faulty within the warranty period, contact CurrentBody.

- CurrentBody reserves the right to repair a defective product.
- CurrentBody reserves the right to replace a defective product with the most compatible product currently available.
- This warranty does not cover damage caused by the following:
 - Accident, misuse, abuse or alteration
 - Servicing by unauthorised persons
 - Connection to incorrect power sources
 - Use with unauthorised accessories

CurrentBody shall not be responsible for any incidental, special or consequential damages resulting from the use of this product. All implied warranties, including but not limited to implied warranties of fitness and merchantability are limited in duration to two years from date of original purchase. This warranty gives you specific legal rights, and you may also have other rights that vary from country to country. Some countries do not allow the exclusion or limitation of incidental or consequential damages, or limitations on how long an implied warranty lasts, so some of the above limitations and exclusions may not apply to you. This warranty is given in lieu of all other warranties, written or oral, whether expressed by affirmation, promise, description, drawing, model or sample. Any and all warranties other than this one, whether expressed or implied, including implied warranties of merchantability and fitness for a particular purpose are hereby disclaimed. If the product is defective within the warranty period please contact CurrentBody.

TROUBLESHOOTING

Problem	Possible Causes	Solution
Light Emitting Indicator turns yellow and flashes continuously.	<ul style="list-style-type: none"> • Alarm for overheating. • There is a problem with the fans. 	<ul style="list-style-type: none"> • The device should be turned off and then powered back up after sufficient heat dissipation. • If the device still flashes after sufficient heat dissipation, please stop using it and contact customer service.
The device cannot be started.	<ul style="list-style-type: none"> • The power supply is not properly connected. • There are some problems with the socket. 	<ul style="list-style-type: none"> • Ensure that the power supply is properly connected. • Change to a different socket.
When placed against the skin, there is no reaction.	<ul style="list-style-type: none"> • There is an overheating problem with the treatment window or the main unit. • There is a gap between the treatment window and the skin. 	<ul style="list-style-type: none"> • After full heat dissipation, adjust the angle of the treatment window until it's flat against the skin. • The skin sensor indicator will be illuminated if the treatment window touches firmly on the skin.
Although the machine can be started, the indicator light does not illuminate.	<ul style="list-style-type: none"> • A failure of the machine has occurred. 	<ul style="list-style-type: none"> • Please discontinue use of this product and contact customer service.
During treatment, the product becomes hot and the main unit becomes warm.	<ul style="list-style-type: none"> • The machine is running overtime. 	<ul style="list-style-type: none"> • The device should be turned off and used after it has completely dissipated.
There is a strange smell during the treatment.	<ul style="list-style-type: none"> • There is a problem with the cleanliness of the treatment window. • The treatment area has not been well shaved, and the remaining hair is burning. 	<ul style="list-style-type: none"> • The treatment window should be cleaned and wiped. • Ensure that the treatment area is shaved thoroughly.
There is a problem with the main unit, light outlet, power adapter, or power cord.	/	<ul style="list-style-type: none"> • Please discontinue use of this product and contact customer service.

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