



GENERAL GUIDELINE



1-2 MONTHS

Within your first three months of using Capillus laser therapy caps, the overall density of your hair may not change significantly. However, you may experience temporary hair shedding. This is normal, as your body is shedding weaker hair allowing newer, stronger hair to grow.

In fact, this shedding is a sign that red-light therapy is working as intended. With the 6-minute daily treatment, you may expect to see noticeable improvements in both the density and quality of your hair.



3-4 MONTHS

After a couple of months of using Capillus laser therapy caps, you may see less shedding and a slight improvement in the density of your hair.



5-6 MONTHS

Around the half-year mark, shedding should have slowed down considerably, and you may see a noticeable difference in both the density and feel of your hair. Ensure optimal growth when you use Capillus' hair care line - specially formulated to help clean and nourish your scalp to promote healthier hair growth.



6+ MONTHS

Regrowing hair with red-light therapy takes time. But, at this point, you should see a difference in your hair and confidence. Integrate Capillus into your daily hair care routine to help maintain the hair you may have grown and all the future hair you could grow. Do not stop here; keep going and growing!

We recommend you take pictures of your scalp before you begin treatment so you can effectively keep track of results.



BEFORE & AFTER

RIVETING RESULTS

No matter your gender, age or hair goals.

FEMALE BEFORE & AFTER



"I'm really happy I found something that worked for me."



"It really changed how I feel about myself."

MALE BEFORE & AFTER



"Capillus is really making wonders in my hair."



"You feel better when your hair looks better."

SETTING UP

Before using the cap, make sure that the battery is fully charged. We recommend leaving it to charge overnight to ensure that the battery is at full capacity.

Once the battery is fully charged, unplug the battery from the charger. Begin your treatment by connecting the cap to the battery. Align the arrows on either cable and gently connect them while making sure they are securely attached. Place the cap on your head and hold down the power button on the battery. You should hear an audible "beep", which indicates that the lasers are on and working.

After 6 minutes, you will hear 3 audible "beeps" from the battery - notifying you that your treatment is complete for the day.

TAKING CARE

To clean your device, disconnect all cables from the laser dome. With an alcohol wipe or damp cloth, wipe and pat dry the inside of the laser dome.

DO NOT APPLY EXCESSIVE PRESSURE TO THE INSIDE OF THE LASER DOME, AS THIS CAN DAMAGE THE DEVICE. DO NOT USE ANY ACETONE OR SOLVENTS AS THIS CAN DAMAGE THE DEVICE AND VOID YOUR WARRANTY.

NOTICE


- Before use, make sure to properly align arrows on the connectors for both the laser dome and battery pack.
- The laser dome will not turn on unless it is placed on the head.
- If used consecutively, allow for 10 minutes before each session to maintain battery health.

FOR MORE HELP ON HOW TO GET STARTED,
VISIT WWW.CAPILLUS.COM/QUICKSTART

BATTERY NOTICE

CHARGE YOUR DEVICE


To ensure optimal usage, charge the battery pack to full (indicated by the 4 LED lights on the battery). We recommend leaving it to charge overnight, as this will give ample time for the battery to charge completely. You can check your battery percentage while undergoing your treatment.

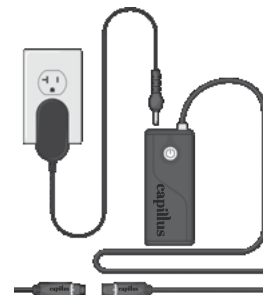

NOTE: YOUR CAP WILL NOT TURN ON WHILE BATTERY IS PLUGGED INTO THE WALL.

CONNECT BATTERY TO CAP

Once the battery is fully charged, gently connect the cable from the battery to the cable from the cap by aligning the arrows by each connector. Make sure the cables are securely connected before turning on the device.

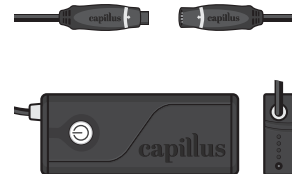

STARTING YOUR TREATMENT

Place the cap on your head. The laser dome will not turn on unless you are wearing the cap.



Then, hold down the power button for a few seconds until you hear an audible beep indicating that you have started your treatment. Treatments last 6 minutes.

AFTER 6 MINUTES...

When your treatment is over, you will hear 3 audible beeps. At that point, the lasers will shut off and you can take off the cap.



DISCONNECT THE CAP FROM THE BATTERY AND STORE IN A COOL, DRY PLACE.

capillus

WARRANTY REGISTRATION

Activate your warranty by visiting www.capillus.com/warranty. Make sure to have your purchase information like price, date of purchase, seller information, and serial number. Your device's serial number can be located on the outside of the packaging and on the rear of the laser device (under the cap). Refer to of your manual for warranty details.

WARRANTY

JOIN THE COMMUNITY
WWW.CAPILLUS.COM

