


WHAT'S IN THE BOX



Curalight (Body)

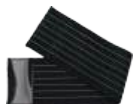

Elastic Strap



USB C Cable

MAINTENANCE AND STORAGE

Cleaning the device

Prior to cleaning, make sure that it is completely unplugged from any power source and is shut off.

Use a damp towel or cloth to wipe away any debris found on the device. **Do not use the damp towel or cloth to clean the LEDs.**

Do not use any liquid or chemical solvents for cleaning, as this may damage the device.

Storing your device.

When not in use, store your device in a safe, dry, and cool place.



BENEFITS OF RED-LIGHT THERAPY



Muscles and joint pains



Enhance blood circulation



Anti-inflammatory



Decrease pain

SPECIFICATIONS

Rating power: 60W

Product size: 400mm x 200mm x 4.95mm

LED quantity: 360 LED diodes (120 x 0.5W)

Net weight: 0.58kg

Wavelength: 630nm & 850nm

Timer: 15-minutes sessions

Power Adapter: 5V / 2.4A (Not included)

Curalight (Body) comes equipped with 630nm, 660nm, and 850nm LEDs. 630nm and 660nm LEDs emit red light, which is visible to the human eye. The 850nm LEDs emit infrared light, which is not visible to the human eye.

WARRANTY INFORMATION

Scan the QR code



CURLPR_23003_R042924

CuraLight BODY

SAFETY INFORMATION

Please read all instructions before using your Curalight device. Keep this user manual for future reference and be sure to follow all warning and safety instructions included.

- Follow all instructions and heed all warnings outlined in this manual.
- Do not stare directly into the LEDs. Prolonged exposure can strain or damage your eyes.
- Use this device only for its intended use and according to the instructions found in this manual. Any other use may result in injury.
- This device is not waterproof. To prevent damage or injury, do not submerge the device in any liquid. If the device does fall into liquid (while the device is operational), unplug from the power outlet.
- DO NOT REACH INTO THE LIQUID TO TURN OFF THE DEVICE.**
- Always plug into a power outlet first, and then plug in the device by using the barrel connector. The device should always be unplugged fully before being plugged into a power outlet.
- DO NOT OPERATE THE DEVICE WHEN IT OR ANY OF ITS COMPONENTS ARE DAMAGED.**
- Do not use the device while asleep.
- Do not use on areas of insensitive skin. Ensure that your skin is dry and free of any chemicals before use.
- This device is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
- If any component is damaged or if any part needs servicing, especially the battery, please contact Curallux.
- Do not lie on top of or place weight onto the device.
- Do not sit on the device or crush the device.
- Do not use the device while operating a motor vehicle.
- This device is not meant to be used by children or pregnant women.
- Do not use any fasteners or tighteners to fasten Curalight to your body.
- Do not wrap device tightly around target area. Allow your skin to breathe during treatment. Ensure that Curalight is fastened securely but not tightly.
- Do not place anything over the device while in use such as a towel, blanket, clothing or any other insulating material.
- This device is not intended to be used on open wounds.
- Do not use device in the hot sun.



GETTING STARTED

Plugging in the device

- Plug-in the **USB C Cable** to a **5V / 2.4A** power adapter (**Not included**).
- Connect **USB C Cable** to the device on the power button.
- When the device is fully charged, the power indicator will turn **red**.



MODES OF OPERATION

- Press **once** and all LEDs will turn on.
- Press **once again** and all LEDs will turn off.



ON



OFF

After you begin your session, the device is auto-programmed to shutoff after 15 minutes.

Curalight (Body) comes equipped with 660nm and 850nm LEDs. 660nm LEDs emit red light, which is visible to the human eye. The 850nm LEDs emit infrared light, which is not visible to the human eye.

AREAS OF USE



- Neck
- Legs



- Shoulders
- Arms



- Upper & Lower Back

Use the device only once or twice a day on the same area. You can use the device on different parts of your body multiple times a day. To ensure the device is used optimally, allow the device to fully cool down before each use.