A man and a woman are smiling and holding a white and blue facial cleansing tool against their cheeks. The man is on the left, and the woman is on the right. They are both wearing white t-shirts. The background is a plain, light color.

**BEAUTYBIO**

**GLOFACIAL**

**HYDRO-INFUSION PORE CLEANSING TOOL**

# A NOTE FROM OUR FOUNDER

“Healthy skin doesn’t have to be complicated. I founded BeautyBio with a commitment to deliver clean, clinical skincare solutions for pro-grade results with zero downtime. Like a vacuum for all of the gunk buried in our pores, our patent pending GLOfacial tool will be your new favorite at-home treatment after the very first, satisfying use - promise. Every time you GLOfacial, turn that routine into a ritual - a minute for you and your skin health. Welcome to the Glommunity!”

Jamie O'Banion  
Founder, BeautyBio

## GOT EVERYTHING?

# GLOFACIAL CONCENTRATE

Clarifying Salicylic and plumping Hyaluronic acid solution helps resurface, hydrate and boost radiance. Niacinamide works to minimize the appearance of congested pores while antioxidant-rich botanical complex with Cica and Green Tea delivers essential amino acids and vitamins for calm, glowing skin.

**INSTANTLY HELPS RESURFACE**

**BRIGHTENS TONE**

**HYDRATES SKIN**

**ALSO SOLD SEPARATELY:**  
0.7 fl. oz. = 8 Facials, 1.7 fl. oz. = 20 Facials

*Unscrew and remove seal before first use.*

# LET'S GET STARTED

**STEP 1**  
Press desired Treatment Tip firmly onto Tool with tab positioned upward.

**STEP 2**  
Slide chamber window and base downward off tool.

**STEP 3**  
Remove Treatment Chamber (left) and pour Concentrate to fill line. Fill remainder with water to dilute.

**STEP 4**  
Reattach Treatment Chamber and slide base back onto tool (note: Waste Chamber does not detach).

**STEP 5**  
Press and hold power button to turn on Blue LED Light mode.

**STEP 6**  
Push button again to start in low suction mode.

**STEP 7**  
Hold skin taut, keep tool upright and glide along damp skin until Treatment Chamber is empty.

**PRO TIPS:** Prime tool on the palm of your hand until concentrate is felt. For best results, we suggest beginning your facial with clean, damp skin. Keep tool moving, do not hover in any one spot.

# THAT'S A WRAP

## STEP 1

Tool will chime when Waste Chamber is full.

## STEP 2

Press and hold button to turn tool off. (Tool will shut off automatically after 10 minutes.)

## STEP 3

Slide chamber window and base downward off tool.

## STEP 4

Open bottom flap of Waste Chamber to empty (note: chamber does not detach).

## STEP 5

Use included brush and warm water to rinse Treatment Tips and chambers (do not submerge tool). Mild soap can be used.

## STEP 6

Dry chambers and treatment tips. Leave Waste Chamber flap open and do not reassemble parts until completely dry.

## STEP 7

Return assembled tool to charging stand.

# CHOOSE YOUR GLOW

## BLUE LIGHT TREATMENT

Press and hold power button to turn on clarifying blue LED treatment mode. Treat with Blue Light for 4 minutes per target area daily.

## PORE-CLEANSING FACIAL

Recommended use 1-2x a week. Treatment is gentle enough to be used on low suction mode daily. Use tool with water alone or combined with GLOfacial concentrate.

Apply gentle pressure to maintain suction and keep tool moving. Do not hover in any one spot or pass over any area more than 1-2 times. Too many passes over an area may cause unnecessary irritation. Some areas of congestion may take several treatments over time to clear. Skin may become temporarily flush after treatment as GLOfacial enhances microcirculation.

*Blue Light Treatment mode does not illuminate chamber.*

# WHAT'S YOUR MODE?

In addition to the blue LED treatment mode, the GLOfacial Tool features 3 suction modes (low/med/high) for customized treatment. We recommend starting with the lowest suction and working your way up. Do not exceed low suction mode on thin or sensitive skin. Do not use directly under eyes or on eyelids.

## BEGIN BY POWERING ON TOOL

(Press and hold button once.)

**LOW**  
Press again for  
low suction.

**MEDIUM**  
Press again for  
medium suction.

**HIGH**  
Press again for  
high suction.

# WHAT'S THAT SOUND?

X1

## BEEPS ONCE

Tool turned on/switching between Blue Light and suction modes.

X3

## 3 BEEPS (WHITE LED FLASHES)

Tool has been tilted too far and is not making optimal contact with skin. Reposition device so it is upright. Do not turn tool upside-down when tool is on.

X3

## 3 TONE CHIME (WHITE LED FLASHES)


Waste Chamber is full and needs to be emptied. Device will not operate again until Waste Chamber is emptied.

# GLIDE PATHS

Follow the natural curves of your face. Remember to hold skin taut, keep moving, and do not hover over one area.

**SUGGESTED GLIDE PATH**

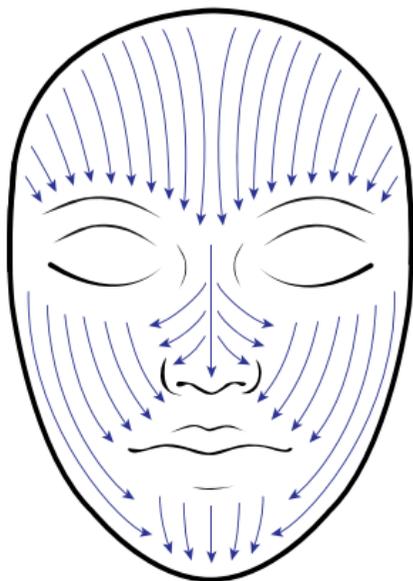**ALTERNATE GLIDE PATH**

**DO** 

- Do make sure your tool is charged.
- Do ensure the tool is clean and dry before use.
- Do choose a treatment mode best suited for your skin.
- Do pull skin taut and glide following facial contours.
- Do hold tool upright and maintain full contact with tip on skin.
- Do empty Waste Chamber immediately after use.
- Do clean and dry chambers after treatment before placing back in charging stand.

**DON'T** 

- Don't overfill Treatment Chamber.
- Don't pour serums or creams into the Treatment Chamber.
- Don't hover in any one spot - keep tool moving.
- Don't pass over any area more than 1-2 times or apply excess pressure with tool.
- Don't hold tool upside down during use.
- Don't use under eyes/on eyelids or shine light directly into eyes.
- Don't let fluid sit in Waste Chamber after use.

# CHARGING

Insert USB charging cable into bottom of charging stand. Connect charging stand to a power outlet (USB power adapter not included). Place tool into charging stand, ensuring gold pins on charging stand are making contact with charging port on back of tool.

Charge GLOfacial tool before first use. Charging time is 2.5 hours. Once fully charged, battery will last approximately 1 month or up to approximately 40 uses. Ensure Treatment and Waste Chambers are empty and tool and charging stand are dry before charging.

## BATTERY INDICATOR LIGHTS

**LOW BATTERY**  
(Red)

**CHARGING**  
(Flashing Green)

**CHARGED**  
(Green)

# FAQ

## 1. What is a GLOfacial treatment?

GLOfacial combines three powerful technologies into one pro-grade treatment at home. Blue LED clarifies skin, gentle facial vacuum technology extracts and exfoliates while hydro-infusion technology flushes pores and infuses skin with nourishing and hydrating ingredients for a deep clean and radiant complexion.

## 2. How often can I use GLOfacial? Do I have to use it at night?

For best results we recommend using GLOfacial 1-2 times a week to maintain results. Start with the lowest suction and working your way up. Some areas of congestion may take several treatments over time to clear. Skin may become temporarily flush after treatment as GLOfacial enhances microcirculation so we recommend treating your skin at night. Blue LED treatment mode can be used for 4 minutes per target area daily.

## 3. Can I use the tool with just water?

You can use GLOfacial with water alone, however you will see best results by treating your skin with the GLOfacial Concentrate.

## 4. Should I wash my face before using the tool?

We recommend beginning with cleansed skin to help prevent makeup, oil or residue from clogging tool and treatment tips.

## 5. Which products can I apply after my treatment to maximize results?

You can follow your GLOfacial treatment with your favorite serums or moisturizers as desired.

## 6. How much pressure should I apply with the GLOfacial tool?

Apply gentle pressure to maintain suction and keep tool moving. Do not hover in any one spot or pass over any area more than 1-2 times. Too many passes over skin may cause unnecessary irritation.

**7. How soon can I start seeing results?**

You will see the debris and oil pulled out of your skin appear in the Waste Chamber immediately after your treatment! Continue using the tool 1-2x a week to maintain your results.

**8. Why are there different tips and do I need to replace them?**

The Large Treatment Tip is ideal for targeting areas like cheeks, forehead and jawline. The Small Treatment Tip targets narrow areas like nose and chin. Over time you may see normal wear and tear. Replace tips as needed.

**9. Is there anywhere on the skin that I should avoid using my GLOfacial?**

We do not recommend using the GLOfacial directly underneath the eyes or on the eyelids. Avoid shining the Blue LED Light directly into your eyes. We also recommend avoiding usage of GLOfacial on any skin irregularities, which include but are not limited to: active acne, open wounds, scars or any raised surfaces including moles. If you have any questions or concerns about specific areas, please contact your physician.

**10. Can I use GLOfacial while pregnant or nursing?**

Always consult a physician before using the GLOfacial tool or GLOfacial Concentrate while pregnant or nursing.

**11. Can I use GLOfacial with fillers or injections?**

We recommend consulting your physician for a personal treatment recommendation. Do not use GLOfacial directly after receiving fillers or injections.

**12. Can I use GLOfacial in the shower?**

GLOfacial may be used in the shower. Avoid submerging tool completely under water. Charging stand is not waterproof or water resistant. Do not use charging stand near water.

**SPECS**

Model Number: GF100

Blue LED: 415nm

Battery: Lithium-Ion Battery 3.7V 1000mAH

Charging Voltage/Current USB: 5V/1.0A

Power Consumption: Approximately 3.7W

Charging Time: Approximately 2.5 Hours

Waterproof Level: IPX7 (Charging stand is not waterproof or water resistant.

Do not use charging stand near water.)

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

**Warning:** Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

**NOTE:** This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules.

These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

## WARNINGS

- Avoid direct eye contact. Avoid use directly under eyes and on eyelids. Do not shine device lights directly into your eyes.
- Do not use on broken skin or skin with severe, nodular or cystic acne.
- Do not use on beauty spots/moles, warts, ingrown hairs, irritated or sunburned skin.
- Avoid passing device over skin multiple times which can cause unnecessary irritation.
- Do not exceed recommended treatment times or usage.
- Do not use with compromised or contaminated liquids.
- Do not attempt to open or repair the device.
- Immediately turn off device if it begins to overheat.
- This tool is not a toy. Keep out of reach of children.
- Do not submerge tool in water. Tool may be used in the shower.
- Do not use it if your skin is light sensitive or if you have a silicone allergy.
- Consult your physician before use if you are pregnant or nursing.
- A medical professional should be consulted before use regarding any condition that requires medical attention.
- Do not ingest.
- Discontinue use if irritation occurs.
- Do not charge tool for more than 24 hours. Remove charging cable from charging stand to continue using as display stand after charging.
- Always use the device in line with the manufacturer's instructions. Improper maintenance, storage or negligence can cause damage and void warranty.

**Precaution:** This device contains a Lithium-Ion Battery (3.7Wh) and may be required to be placed in carry-on baggage for travel. Check your local airline regulations for guidance on how to travel with your tool. Do not dispose of device as household waste.

If you have any additional questions, please contact a customer care specialist at **877-945-5274** or **info@beautybio.com**. For additional information, please visit **BeautyBio.com**

## RECYCLING OF ELECTRONIC PRODUCTS AND BATTERIES

Electrical and electronic equipment (EEE) contains materials, components and substances that may be hazardous and present a risk to human health and the environment when waste and electronic equipment (WEEE) is not handled correctly.

Equipment marked with the below crossed-out wheeled bin is EEE. The crossed-out wheeled bin symbol indicates that WEEE should not be discarded together with unseparated household waste, but must be collected separately.

Users of EEE must not discard WEEE together with household waste. Users must follow local recycling regulations to reduce adverse environmental impacts in connection with disposal of WEEE and to increase opportunities for reuse, recycling and recovery of WEEE.

Discarded household electrical and electronic equipment and batteries can be returned to consumer product collection points free of charge. The nearest collection point can be found easily from this website: [www.recyclenow.com](http://www.recyclenow.com).

In addition to this recycling drop off locations BeautyBio is a member of the Distributor Take Back Scheme (DTS). For more information on recycling options in the UK please go to: [www.recycle-more.co.uk](http://www.recycle-more.co.uk)

Where possible, used products should have batteries removed to allow for those batteries to be recycled.

Beauty Biosciences, LLC is registered as a WEEE producer with the UK Environment Agency (WEEE registration number: **WEE/AX5638XC**).

Packaging is crafted with reduced and recyclable materials. Thank you for supporting a sustainable future.

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