

osmo

.....

EYEREST USER MANUAL





[RATED VOLTAGE] DC 5V
[BATTERY CAPACITY] 1000MAH
[PRODUCT WEIGHT] 376G
***1 USB DATA CABLE*1**

Product Introduction

The intelligent eye massager incorporates cutting-edge microcomputer chip control technology to perform kneading, vibration, and heat compression on the eyes. This multi-faceted approach is designed to alleviate eye muscle tension, delivering an all-encompassing sense of comfort and relaxation.


This product caters to individuals who extensively utilize computers, mobile phones, and other electronic devices. To begin using it, please follow these steps:

1. Insert the USB plug of the USB data cable into the USB socket of the power adapter (5V-1A).
2. Connect the power adapter (5V-1A) to a power socket.
3. Finally, insert the Micro USB plug of the USB data cable into the Micro USB socket of the product.

A friendly reminder: While charging, the indicator light will flash. Once fully charged, the indicator light will remain illuminated. During the charging process, the product cannot be activated.

After charging is completed, kindly disconnect the power adapter. If the product is charging while in use, it will automatically cease operation to ensure your safety.

Product Structure

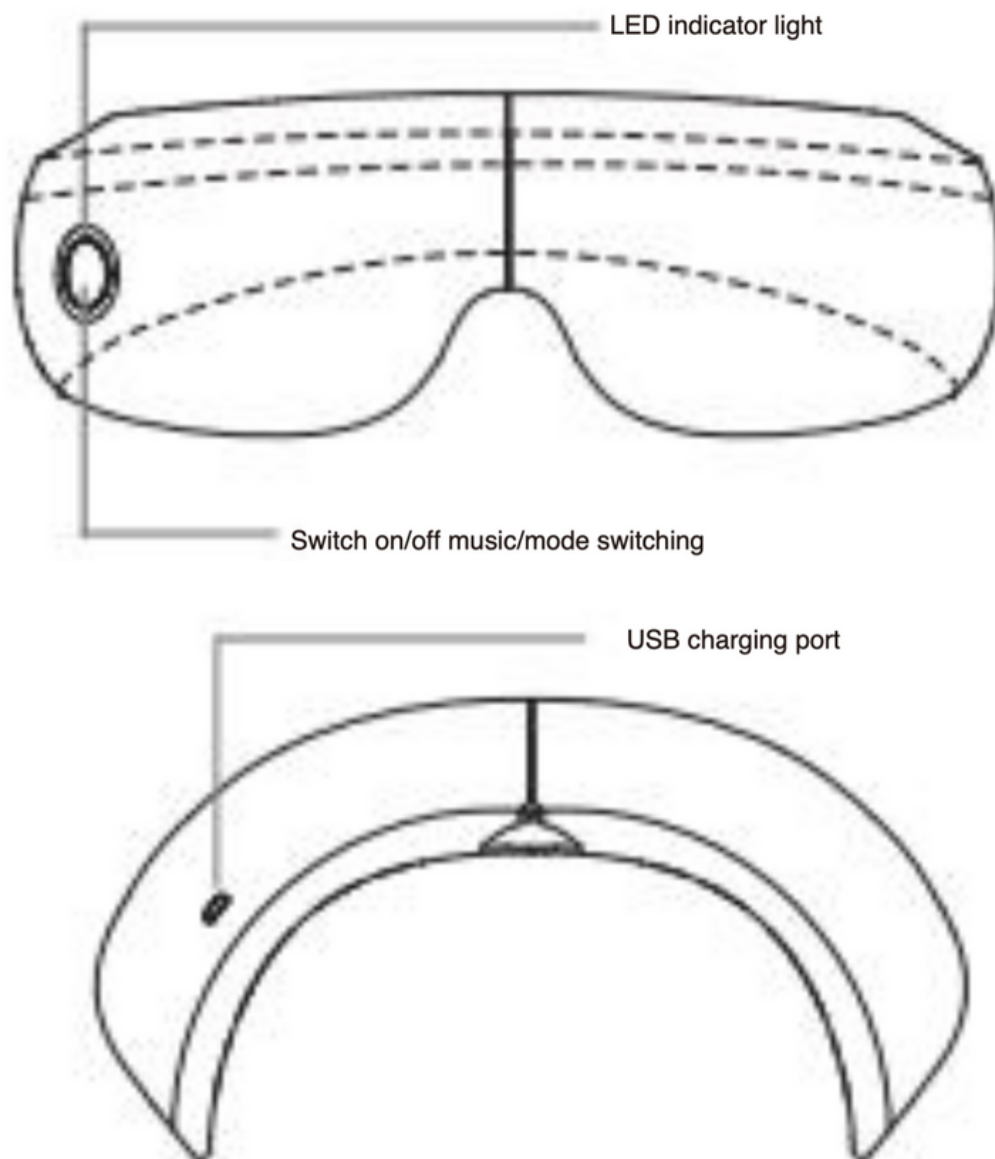

Product Structure

The product includes:

LED indicator light

On/Off switch for music and mode switching

USB charging port

Instructions

Power On/Off

- To activate the product, press and hold the "Power" button for 2 seconds. The indicator light will illuminate, and the machine will begin operation. The default mode is automatic, and it will run for 15 minutes while playing music.
- While the product is in operation, press and hold the "Power" button for 2 seconds to turn it off. Alternatively, the machine will automatically shut down after 15 minutes of use.

Switching Modes

- For mode selection, briefly press the "Mode" button. A single press switches to the second mode, while a double press switches to the third mode, and so on.

Music Control

- To manage the music, press the "Music" button twice continuously in under one second. This action will either turn off or turn on the music.

Mode Introduction

- Automatic Mode: Combines air pressure massage, hot compress massage, vibration massage, and music.
- Eye Protection Mode: Integrates air pressure massage, hot compress massage, and music.
- Vitality Mode: Features air pressure massage, vibration massage, and music.
- Comfort Mode: Focuses on hot compress massage and music.
- Sleep Mode: Involves air pressure massage and music.

Instructions

Bluetooth Connectivity

- Once the eye protection device is turned on, Bluetooth is activated by default. To establish a Bluetooth connection, please enable Bluetooth on your mobile phone and search for devices. The default Bluetooth name is "Eye Protector." Upon successful connection, you can play music from your phone. After disconnecting Bluetooth, the device will play built-in meditation music.

Transportation and Storage

Storage Environment

- Store the product in a well-ventilated, dry room with a normal temperature and no corrosive gases.

Transportation

- This product is equipped with basic shock-proof features, making it suitable for transportation via aviation, railway, highway, or ship. Ensure it is safeguarded against rain, exposure to the elements, inversion, collisions, and other adverse conditions.

Cleaning and Storage

- When cleaning the product, use a dry, soft cloth, and remember to turn off the machine. Avoid using corrosive cleaning agents such as benzene or thinner, as they can lead to issues like fading, surface aging, and structural damage.
- For long-term storage, place the product in a well-ventilated, dry room, and consider using a protective cover. If the product will be stored for over three months, it is recommended to unpack and check it, charging it for at least half an hour.

Safety Precautions

- 1. Eye Conditions:** Refrain from using this product if you exhibit symptoms of eye disease or have undergone eye surgeries. Individuals suffering from major eye ailments like glaucoma, cataracts, or retinal detachment should avoid its use.
- 2. Cleaning Instructions:** Do not attempt to clean this product directly with water. Instead, use a cotton cloth soaked in either alcohol or water for cleaning.
- 3. Heat Sensitivity:** This product features a heating surface. Exercise caution when using it on individuals who may not be sensitive to heat, such as minors, people with disabilities, and the elderly. In such cases, adult supervision is necessary.
- 4. Sensitive Skin:** For those with sensitive skin, allergic dermatitis on the face, or individuals who have plastic materials on their face, it is advisable to use this product under the guidance of a medical professional.
- 5. No Disassembly:** Do not attempt to disassemble this product on your own. It contains batteries, and placing it in high-temperature environments or near open flames may lead to hazardous situations, including the risk of explosion. Always follow proper usage guidelines.
- 6. Blocked Line of Sight:** Using this product will obstruct your line of sight. Avoid using it while engaged in activities such as walking, driving, and other tasks that require unobstructed vision.
- 7. Post-Use Blurred Vision:** After using this product, it is normal to experience a brief period of blurred vision as your eyes transition from a relaxed to a tense state. To alleviate this, close your eyes for 1-3 minutes, and your vision will gradually clear.

Troubleshooting

Issue:

- The device automatically shuts down after being turned on or cannot be turned on.
- The charging indicator light does not illuminate.

Possible Causes:

Low battery charge.

Solution:

Please charge the device.

Product Features

This product features automatic timing, far-infrared heating, voice broadcasting, portability, rechargeable battery, Bluetooth music connectivity, multi-frequency vibration, and 180-degree folding capability.

Important Preparations and Charging Before Use

Preparation Before Use:

Align the eye protection device with your eyes, wear the device, adjust the elastic band to a comfortable size, and position the device up, down, left, and right for the most comfortable fit.

Prior to use, remove glasses, contact lenses, hairpins, hair accessories, and the like to prevent discomfort. Avoid using the product with wet hands, and take care not to drop or step on it.

Risk Considerations:

This product serves as a personal health massage device and is not a substitute for medicines or medical equipment.