

osmo

TENS UNIT USER MANUAL



This product is not a medical device product.

Safety precautions:

Do not use with the following medical electronic devices:

1. Use of pace makers and implantable medical electronic devices;
2. Use of artificial cardiopulmonary pacemaker and life extension medical electronic devices;
3. Using ECG devices and wearing medical electronic devices;
4. The placement of the massage patch should be avoided from the heart and chest of the human body. No electrical pulse to affect the heart causing discomfort;
5. Skin abrasions and suppurative Inflammation, use with caution; When getting a massage from one of the above, Massage patch to avoid the wound, so as not to sting and aggravate the wound;
6. It is recommended to be used under the guidance of the seller. Patients, children and children without Independent ability are recommended. Patients with heart disease, diabetes and liver cirrhosis should not be allowed.
7. Patients with implanted metal devices are prohibited to use; If the above medical electronic devices are used at the same time, it may be dangerous due to wrong operation.

Use with caution where:

1. Malignant tumor;
2. Infectious diseases;
3. Menstruating and pregnant women;
4. A heart attack should not be placed near the heart.
5. Fever above 38 degrees;
6. Abnormal skin sensitivity;
7. In the following parts of the body: mouth, eyes, perineum, skin trauma, Do not use this product;
8. Do not use it together with other massage machines or similar Instruments or application of medicines.
9. Do not use this product on an empty stomach;
10. Please be careful when you feel uncomfortable or dizzy. This product belongs to massage Instrument, not to medical device products.

Functional Features

- 1) Apply digital technology to transform physical factors needed by human body into precision which is then massaged by a computer according to the needs of the human body, Ensure the accuracy of massage, greatly improve the efficacy of manage.
- 2) With a large Screen LED, covering strength, function, mode, massage. Time. Chinese and English display system, meridian diagram of human body.
- 3) Massage patch is used to cover a larger area.
- 4.) Various Massage Functions. Increase and coordinate body resistance, effectively slow Body aches and pains caused by various external factors.
- 5) wonderful wave form combination, so that you can truly experience the simulation of acupuncture, pounding The real feeling of stroke, massage, cupping, massage and scraping
- 6) LED meridian diagram for massage. According to their own needs, and by the meridian diagram Guide, Select the sore and uncomfortable parts for massage. easy to operate.

Product Features

This product USES the high-tech computer digital technology, The imitation hand like massage Movement, in integration of advanced ergonomic theory, similar to the human hand skin surface press Rub patch, thin 2 pieces stick on the place needing manage, accelerate blood circulation, fast quickly eliminate fatigue, neck pain, and stiffness, and quickly relieve the pressure or work. Bring back the good cheer.

This product models exquisite, compact. light., can be used in a warmly family. Yes, it can also be used in cars and offices, once to three times a day, per day Five to fifteen minutes all types of muscle soreness, back soreness. and shoulder pain, Neck pain, lag pain. etc., simulated hand massage, fatigue muscles quickly. Rejuvenate. Bring you an unexpected sense of ease and make your work come true.

Instructions:

1. Open the battery cover, install the positive and negative electrodes of the battery and close the battery cover.
2. Insert the button line into the output hole at the top of the product, and then connect the button line with the patch. Tear off the surface film of the patch, and then stick the massage patch to the place to be massaged, note that two massage stick to need to be placed on a person at the same time.
3. Start using the product.
4. This product has two power supply modes: power adapter and battery. You need to take the battery out.
5. Unplug the power source when not in use .
6. USB cable can also be directly connected to the computer for use.


Functions:

1. Button line socket: connect the massage patch
2. Liquid crystal display screen: Intuitive synchronous display of the machine state.
3. Power on/key: power on/power off function.
4. Display key in Chinese/English: Press this button to select LCD display in English
5. Function selection key (mode key): this machine has 8 massage functions: Acupuncture, beat, massage, cupping, scraping, slimming, immune regulation, according to this: Click once to change a massage mode, 8 modes can be selected successively, open the machine automatically enters the needle mode:


6. Timing Key: each time you press this key, you can increase timing for 5 minutes, and timing for 15 minutes at most, Automatic timing for 15 minutes after starting:
7. Intensity Increase Key: Increase the intensity of massage, Each time you press the key, it will increase intensity until it reached the maximum of 15
8. Intensity for the weakening key: weaken the massage intensity, each button once reduce a strong degrees until 0.



Refer to the schematic diagram of massage points

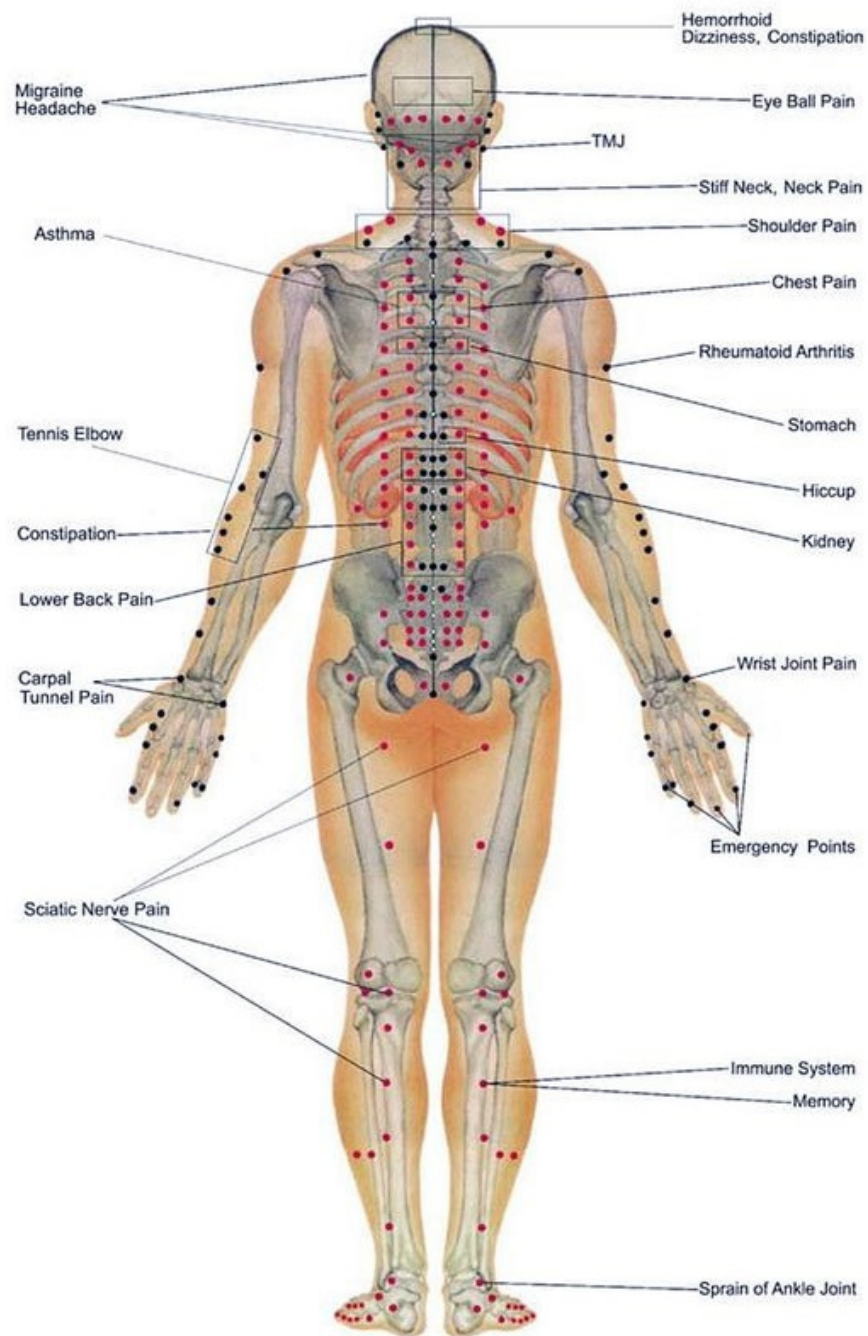

Refer to the schematic diagram of massage points

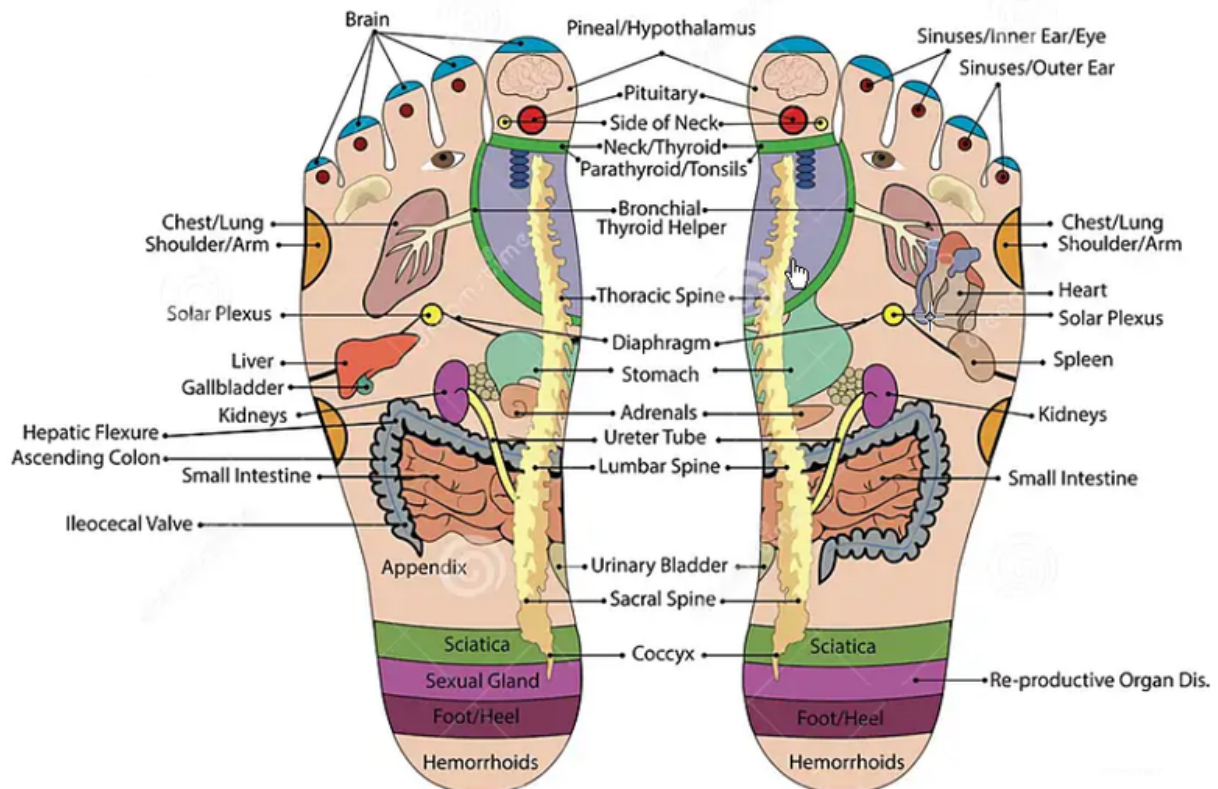

Preparation Before Using the Device

Load Cell

1. Gently open the battery cover
2. Load 3 selection of no.7 battery, and pay attention to the positions of positive and negative electrodes.
3. Close the lid of the box

Prepare Massage Patch

1. Clean the area that needs to be massaged to remove oil stains and on the skin surface.
2. Gently tear off the protective paper on the surface of the massage patch; Protect paper from the abandonment, massage T place the tape on the protective paper after use

Massage Position

1. Dry the area to be massaged with a wet towel to remove the grease and statins Apply the massage patch to the part of requiring massage
2. 2 pairs of the massage patches should be attached to the massage area at the same time, otherwise it cannot be carried out massage patch does not need to placed symmetrically.

After Sales and Warranty:

1. If any performance failure occurs 14 days upon receipt consult the customer service team to check the coverage.
2. 2 year warranty replacement for any damages or faulty items; Damaged due to improper use and not following the storage method is not included

Product Maintenance and Storage Method

1. Clean the main engine; wipe it with a soft cloth, do not wet the main engine, do not use wet wool cloth, Wipe the main engine to avoid the water infiltration, do not use diluents, gasoline or other easy to develop liquids.
2. Clean and massage the patch; the dirt on the patch can be wiped with a wet cloth or washed with water/ Wash immediately and dry to be used.
3. Storage: Place the massage machine out of reach of the children, Keep it away from humidity, high temperature, and direct Sunlight.

Transportation:

The product is equipped with simple anti-vibration facilities, suitable for air, railway, highway and ship transport. Rain and snow splashing, inversion, and collision should be avoided.

Storage:

- If the storage period exceeds 6 months, the box should be unpacked for inspection and electrified for half an hour; not close to the ground, walls and roof
- Indoor ventilation should be good, avoid strong sunlight and corrosive gas erosion

Product Specification

Power Supply: 3 DC 4.5V Dry Batteries

Frequency: 1--330 Hz

Power: 17m A (resistance load = 1K)

Host Size: 1.5.3cm*7.0cm*2.8cm

Weight: 90g

Troubleshooting Guide

Failure State	Why	To deal with
Load the Battery LCD does not display	<p>Whether there is any foreign body in the battery box</p> <p>Whether the battery is exhausted or installed upside down</p> <p>Whether the battery is in the bad contact with the battery</p>	<p>Check and Clear</p> <p>Replace or place the battery correctly</p> <p>Use the right tool to pick out the pole</p>
Load battery liquid crystal flash but does not flash	It is not malfunction	Remove the battery after 4 seconds
The liquid crystal is normal, no sensation	<p>The plug is not correctly inserted into the hole</p> <p>Whether the massage patch is in normal contact with the skin correctly</p>	<p>Correct Plug</p> <p>Button the thread and massage the patch</p>
The stimulation becomes weak or has a numbing sensation	<p>whether the massage patch is in the contact with the skin correctly</p> <p>Massage patch surface is dirty</p> <p>Whether the body surface has grease</p>	<p>Place the massage patch firmly on the skin</p> <p>Clean or replace the massage patches</p> <p>Clean the surface</p>
In massage, the output is stopped before 15 minutes are used	Whether the battery is exhausted	Replace the battery
The skin becomes red, itchy or irritated at the massage site	The skin is not sensitive to sticky pads	Check the skin for any history of allergy. For minor allergies, only change the massage position or shorten the massage time. For several allergies, stop using it